


58 Litre Toploading Autoclave

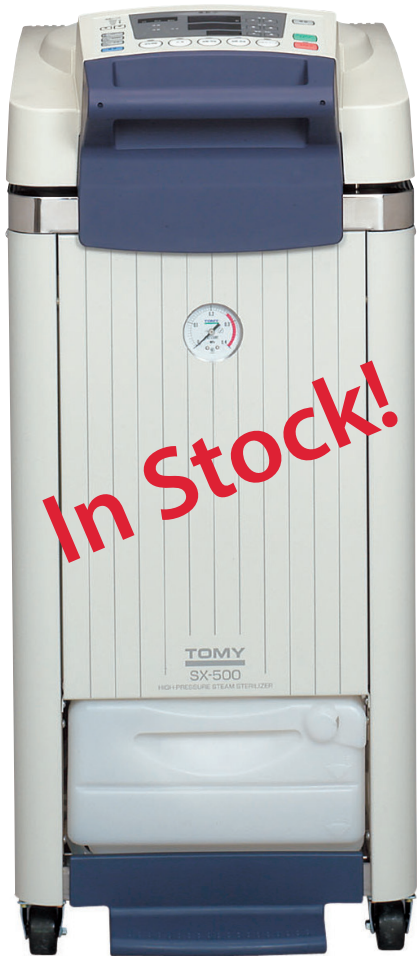

Only \$8,675

Usually \$14,726

Code: TOMYSX-500E
Applies to current stock only!

Supplied with two baskets &
2-year warranty

Having castors for mobility and an integrated condensate bottle, the autoclave requires no plumbing so can be moved to where it is needed!

- Three operating temperature ranges: Warming 45 to 95°C, Heating 45 to 104°C and Sterilizing 105 to 135°C
- Grade 304 stainless steel chamber for easy cleaning and durability
- Variable speed exhaust and cooling fan to reduce run time.
- Basket Dimensions: 300mm Ø x 181mm H
- Five sterilization profiles to select from (in addition to normal sterilization)
- Safety devices include: water level sensor, current leakage breaker, lid interlock, overheat prevention, over pressure prevention, open temperature sensor detection & safety valve.
- Digital temperature & time display with analogue gauge for pressure.
- Integrated (removable) condensate bottle to trap condensate from chamber during venting.



Simple lid opening using
either hand or foot.

Good range of accessories including punched metal stainless steel baskets, solid stainless steel baskets, hybrid basket: solid lower third with punched metal stainless steel upper, Durham tube stand, printer, external sensor (wandering probe) and preventative maintenance contracts.